



DEEP2

SYNERGETIC SKIN REJUVENATION

MULTI-FUNCTIONAL FACIAL PLATFORM

The deep2 dermo-facial platform provides three technologies that, combined with each other and with the met2morph cosmeceutical line, work synergistically to encompass the complete range of modern, professional facial treatments.

The deep2 technology incorporates hydrodermabrasion, draining / remodelling, cupping massage and bio-stimulating LED therapy, which respond to the needs and characteristics of each type of skin.

PURIFIES, NOURISHES & BALANCES

Removes impurities and dead cells, freeing the pores and reducing their appearance. Accelerates the entry of cosmeceutical complexes that nourish, balance and activate skin metabolism.

DECONGESTS & PROTECTS

Stimulates the elimination of retained liquids, freeing the skin of toxins and unwanted metabolic waste. Restores a supple, luminous skin protected by antioxidants, nutrient peptides and trace minerals that rebuild the barrier function.






GLIT : HYDRODERMABRASION

The GLIT handpiece uses diamond tip exfoliation that allows a deep cleansing of the skin, performing a gentle exfoliation that removes impurities from the pores and prepares the epidermis for the infusion of the active ingredients of the treatment.

Its suction-infusion technology allows the application of met2morph cosmeceutical complexes simultaneously with the removal of the exfoliation residues, leaving the skin clean, luminous and ready for the next steps of the deep2 treatment.

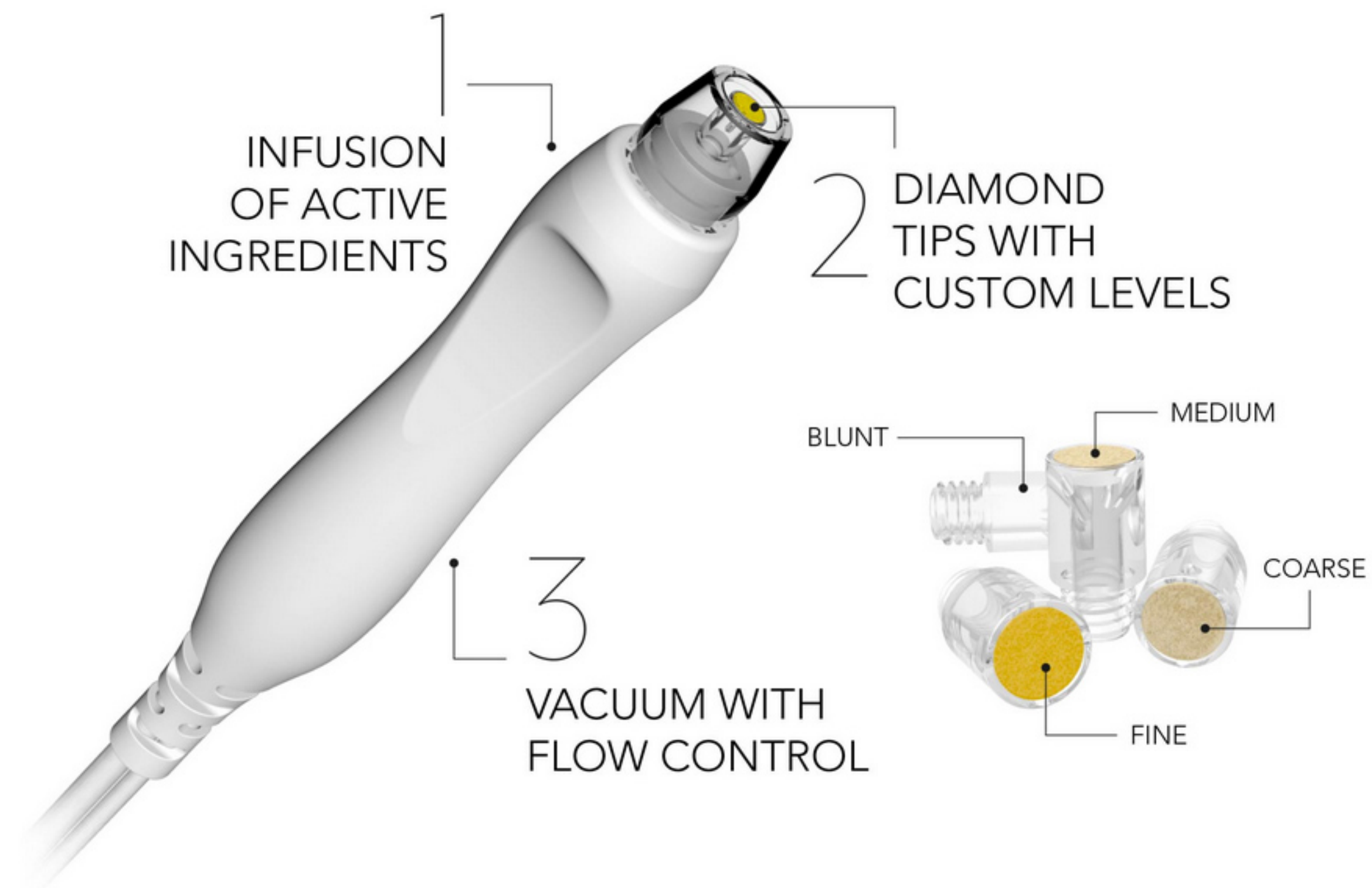

CALYX : CUPPING MASSAGE

With its precise design, the CALYX applicator allows you to perform a relaxing cupping massage that reshapes, drains and decongests the skin according to its specific needs.

- Stimulates blood circulation, improving the supply of oxygen and nutrients to the skin.
- Acts on the lymphatic circulation triggering drainage of the tissue and detoxification.
- Used as a remodelling treatment, promoting the generation of collagen and elastin.






MODES OF USE



CONTINUOUS
REMODELLING
RESULTS



PULSATING
DRAINING &
CIRCULATION

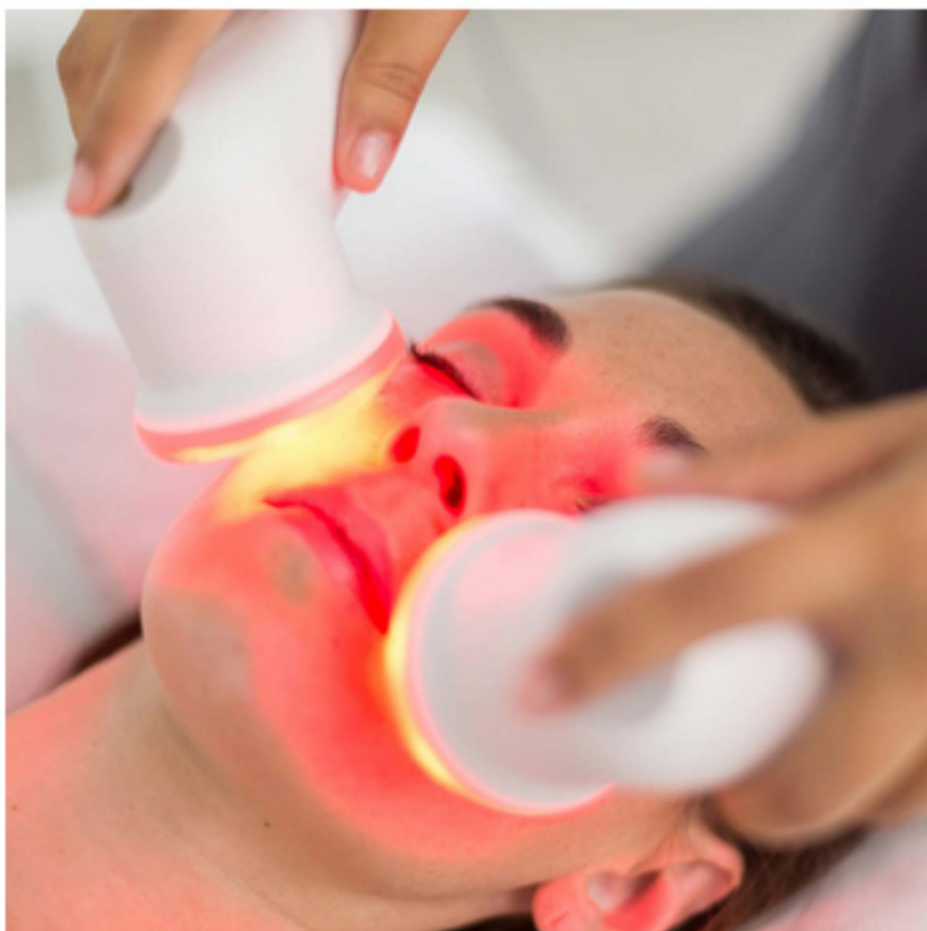

LUX : LED THERAPY

Two simultaneous applicators of bio-stimulant LED therapy that, used together with the CRYBALID mask, exert a calming effect at the end of the deep2 treatment. LUX chromotherapy, with red, green and blue light modes, activates the skin's fundamental repair processes.

LUX : LED BENEFITS

- Calming effect / Balances sebum
- Modulates the microbiota
- Activates repair processes
- Stimulates the synthesis of dermal matrix

LED COLUR SPECTRUM

● 468NM : BLUE LIGHT

Inhibits bacterial proliferation and improves the microenvironment of the skin.

● 630NM : RED LIGHT

Stimulates the production of collagen and elastin, improving skin hydration.

● 468NM : GREEN LIGHT

Unifies skin tone for a more luminous appearance. Stimulates the synthesis of the dermal extracellular matrix.



PROFESSIONAL COSMECEUTICALS

Science, nature and beauty are the elements that inspire met2morph, the cosmeceutical line of the deep2 universe that was created to enhance the health and balance of the skin.

We believe in healthy skin, which transforms and changes, revealing itself to be beautiful at every moment and at every age. Giving us security and confidence, reflecting our essence. For this reason, its balance is important and is unique for each skin condition.



UNIQUE FORMULATION

Expert cosmetologists and facial specialists have created unique formulations with more than 75 state-of-the-art active ingredients that infuse the skin, nourishing it, balancing it and protecting it.

met2morph was created as a perfect symbiosis between cosmeceutical science and deep2 technology. Together, they create a universe of treatment experiences that allow total customisation, responding to the changing needs of the skin.

Because all skins are different and their needs change, deep2 & met2morph allow the total customisation of the treatment and its adaptation to the skin's evolution.



SEE THE RESULTS



SPECIFICATIONS

Supply voltage	230 V AC 50 Hz 115 V AC 60 Hz (optional)
Maximum consumption	150 VA
Interface	5.7" colour touch screen
IP protection rating	IPX0
Classification	Clase I type B
Fuse protection	2 x T1A / 250VAC
Dimensions	18 cm (w) x 34 cm (l) x 40 cm (h)
Weight	8.5 Kg
HYDRODERMABRASION FEATURES:	
Applicator tip diameter	6mm
Applicator tip thicknesses	0, 30, 45 and 60 microns
Maximum pressure	20 cmHg
CUPPING CHARACTERISTICS:	
Suction level	1-10 cmHg
LED THERAPY FEATURES:	
RED wavelength	630nm
BLUE wavelength	468nm
GREEN wavelength	523nm